

## ✓ BPL\_YEAST\_AIR\_Fedbatch script with PyFMI

The key library PyFMI is installed.

After the installation a small application BPL\_YEAST\_AIR\_Fedbatch is loaded and run. You can continue with this example if you like.

```
!lsb_release -a # Actual VM Ubuntu version used by Google
```

```
No LSB modules are available.
Distributor ID: Ubuntu
Description:   Ubuntu 22.04.3 LTS
Release:      22.04
Codename:     jammy
```

```
%env PYTHONPATH=
```

```
env: PYTHONPATH=
```

```
!wget https://repo.anaconda.com/miniconda/Miniconda3-py310_23.1.0-1-Linux-x86_64.
!chmod +x Miniconda3-py310_23.1.0-1-Linux-x86_64.sh
!bash ./Miniconda3-py310_23.1.0-1-Linux-x86_64.sh -b -f -p /usr/local
import sys
sys.path.append('/usr/local/lib/python3.10/site-packages/')
```

```
--2024-01-18 10:10:42-- https://repo.anaconda.com/miniconda/Miniconda3-py310
Resolving repo.anaconda.com (repo.anaconda.com)... 104.16.130.3, 104.16.131.3
Connecting to repo.anaconda.com (repo.anaconda.com)|104.16.130.3|:443... conn
HTTP request sent, awaiting response... 200 OK
Length: 74403966 (71M) [application/x-sh]
Saving to: 'Miniconda3-py310_23.1.0-1-Linux-x86_64.sh'
```

```
Miniconda3-py310_23 100%[=====>] 70.96M 162MB/s in 0.4s
```

```
2024-01-18 10:10:42 (162 MB/s) - 'Miniconda3-py310_23.1.0-1-Linux-x86_64.sh' :
```

```
PREFIX=/usr/local
Unpacking payload ...
```

```
Installing base environment...
```

Downloading and Extracting Packages

Downloading and Extracting Packages

```
Preparing transaction: done
Executing transaction: done
installation finished.
```

```
!conda update -n base -c defaults conda --yes
```

```
!conda --version  
!python --version
```

```
conda 23.11.0  
Python 3.10.13
```

```
!conda install -c conda-forge pyfmi --yes # Install the key package
```

```
Preparing transaction: done  
Verifying transaction: done  
Executing transaction: done
```

## ✓ Notes of BPL\_YEAST\_AIR\_Fedbatch

Now specific installation and the run simulations. Start with connecting to Github. Then upload the two files:

- FMU - BPL\_YEAST\_AIR\_Fedbatch\_linux\_jm\_cs.fmu
- Setup-file - BPL\_YEAST\_AIR\_Fedbatch\_explore

```
%%bash  
git clone https://github.com/janpeter19/BPL_YEAST_AIR_Fedbatch  
  
Cloning into 'BPL_YEAST_AIR_Fedbatch'...
```

```
%cd BPL_YEAST_AIR_Fedbatch  
  
/content/BPL_YEAST_AIR_Fedbatch
```

## ✓ BPL\_YEAST\_AIR\_Fedbatch - demo

Author: Jan Peter Axelsson

This notebook demonstrate yeast fedbatch cultivation. We look at impact of changes in the glucose feeding. We also take a look at tuning of the DO-control system. Both liquid- and gasphase are included in the model. The culture growth and metabolism are formulated in relation to the respiratory capacity [1] and the model is expanded to describe also the gas phase as well as the culture heat production [2]. The model was derived mainly from continuous culture data but proved to capture dynamic aspects well of ethanol production and consumption [3].

Interaction with the compiled model as FMU is mainly through the simplified commands: `par()`, `init()`, `newplot()`, `simu()` etc. The last simulation is always available in the workspace and called 'sim\_res'. The command `describe()` brings mainly up description information from the actual Modelica code from the FMU but is complemented with information given in the dedicated Python setup-file.

The idea is to demonstrate how simulations and varying conditions can provide some process insight that can support the experimental work. I hope that at the end of this session you are ready to formulate your own questions you want to address with simulations - and you can just go on in this notebook! Just press the field "+Code" in the upper left part of notebook interface and you get a new "cell" where you write your own code. You can copy and paste from cells above using ctrl-c and ctrl-p as usual and edit the cell. When you are ready to execute the cell just press the "play button" to the left in the cell or press shift-enter as in "ordinary" Jupyter notebooks.

After a session you may want to save your own notebook. That you can do on your Google Drive account and I refer to Colab instructions for how to do this. It is easy.

Enjoy!

```
run -i BPL_YEAST_AIR_Fedbatch_D0control_explore.py
```

```
Linux - run FMU pre-compiled OpenModelica 1.21.0
```

```
Model for bioreactor has been setup. Key commands:
```

- par()            - change of parameters and initial values
- init()          - change initial values only
- simu()          - simulate and plot
- newplot()      - make a new plot
- show()          - show plot from previous simulation
- disp()          - display parameters and initial values from the last simulation
- describe()     - describe culture, broth, parameters, variables with values/units

```
Note that both disp() and describe() takes values from the last simulation
```

```
Brief information about a command by help(), eg help(simu)
```

```
Key system information is listed with the command system_info()
```

```
%matplotlib inline
plt.rcParams['figure.figsize'] = [36/2.54, 30/2.54]
```

## ✓ About the process model

We can get information about the process, liquid- and gas-phase by the command describe(). This command can also be used to bring up information about a specific variable or parameter. However, you should use describe() after a simulation to get the values used during the simulation.

```
describe('culture'); print(); #describe('liquidphase'); print(); describe('gaspha
```

```
Saccharomyces cerevisiae - default parameters for strain H1022
```

The model of the process has parameters both for culture, gas\_liquid\_transfer, as well as feeding procedure. The parameters that are available for changes you find by the command `disp()` and you get a long list and you change by them by command `par()`. The model has even more parameters in the background but not made available for interaction.

## ✓ First simulations - adjusting start of substrate feeding


```
# Culture parameters and others at default values
par(q02lim=0.0069)

# Process initial conditions
init(V_0=4.5, VG_0=4.5*30, VX_0=4.5*2, VE_0=4.5*0)

# Feed profile
par(t_start=9, F_start=0.044, mu_feed=0.20, F_max=0.18)

# DO-control parameters
par(K=10, Ti=0.5)

# Simulate and plot
newplot(title='Yeast fedbatch cultivation', plotType='Overview')
simu(20)
```

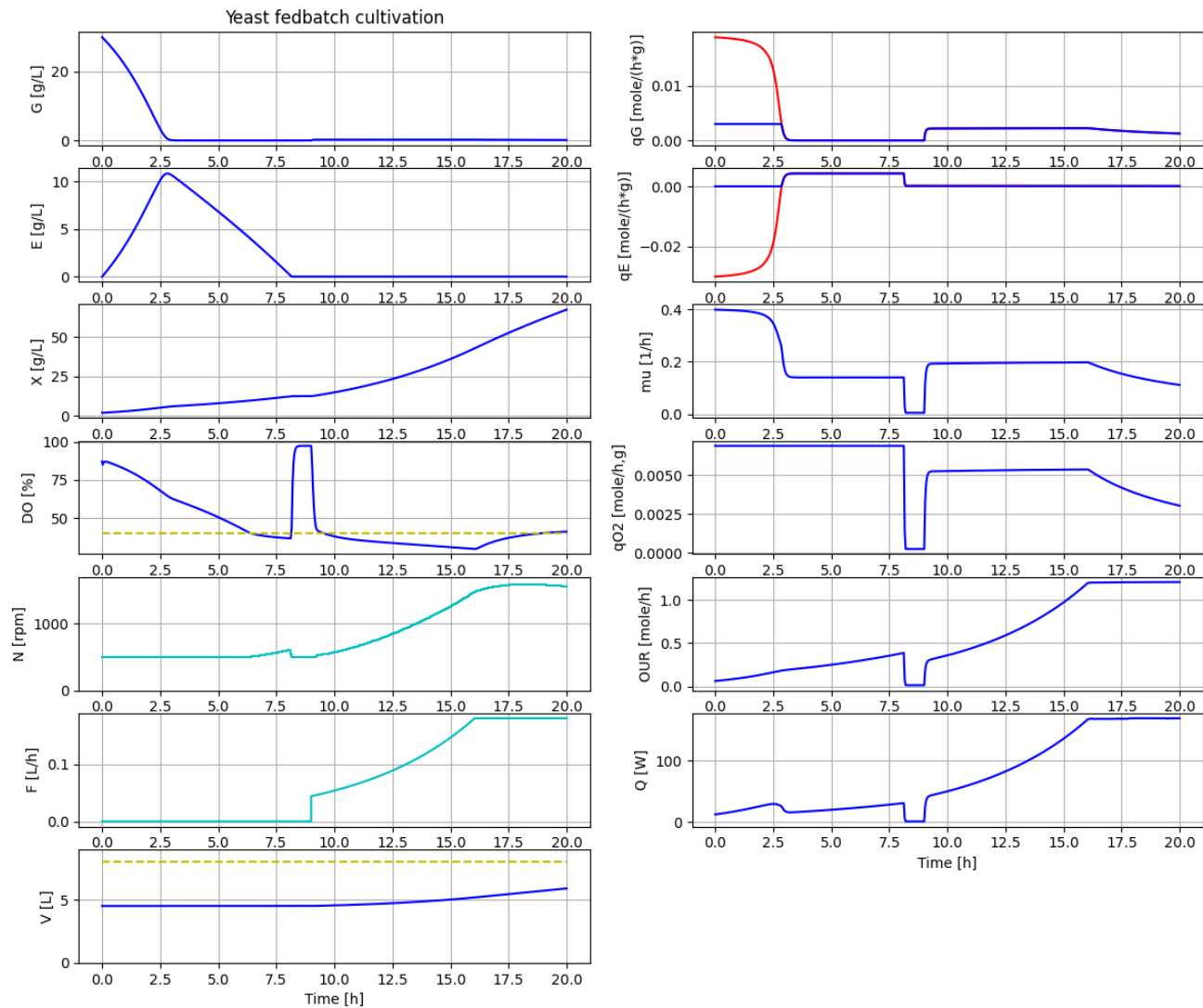

Now we can get value of broth volume as well as the headspace and values are the last ones in the simulation

```
describe('bioreactor.V')
```

```
Reactor broth volume : 5.892 [ L ]
```

```
describe('bioreactor.V_gasphase')
```

```
Volume of the gas phase : 2.108 [ L ]
```


```
# Take a look at the parameters available to adjust the dosage scheme  
disp('dosage', decimals=4)
```

```
mu_feed : 0.2  
F_0 : 0.0  
t_start : 9.0  
F_start : 0.044  
F_max : 0.18
```

```
# Let us start the feeding just after the batch phase has ended and keep other pa  
par(t_start=8.1)
```

```
# Simulate and plot  
newplot(title='Yeast fedbatch cultivation', plotType='Overview')  
simu(20)
```



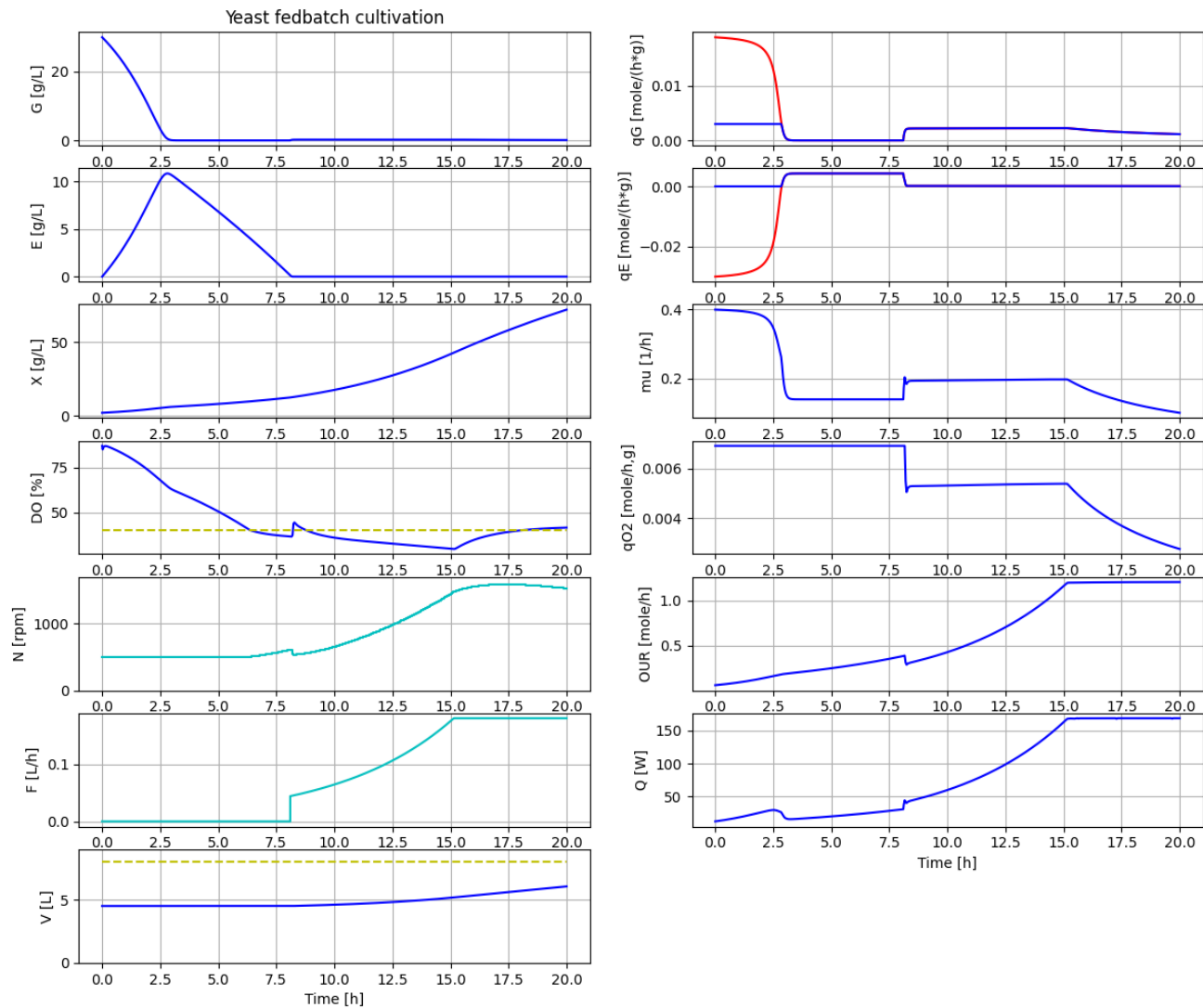


The increase of DO to about 50 % at end of batch phase should be possible to detect easily. This simulation is more realistic and we use these settings from now on.

## ✓ DO-control - tuning of PI-regulator parameters

Let us focus on the DO-control system and choose a more limited plotType. We study the impact of PI control parameters and see if we can decrease the control error without losing stability.

```
# Let us take a closer look at the DO-control system and try to make control error  
newplot(title='Yeast fedbatch cultivation - DO-control - increase K', plotType='F'  
for value in [10, 20, 40]: par(K=value); simu(20)
```

```
# Reset K to the original value  
par(K=10)
```

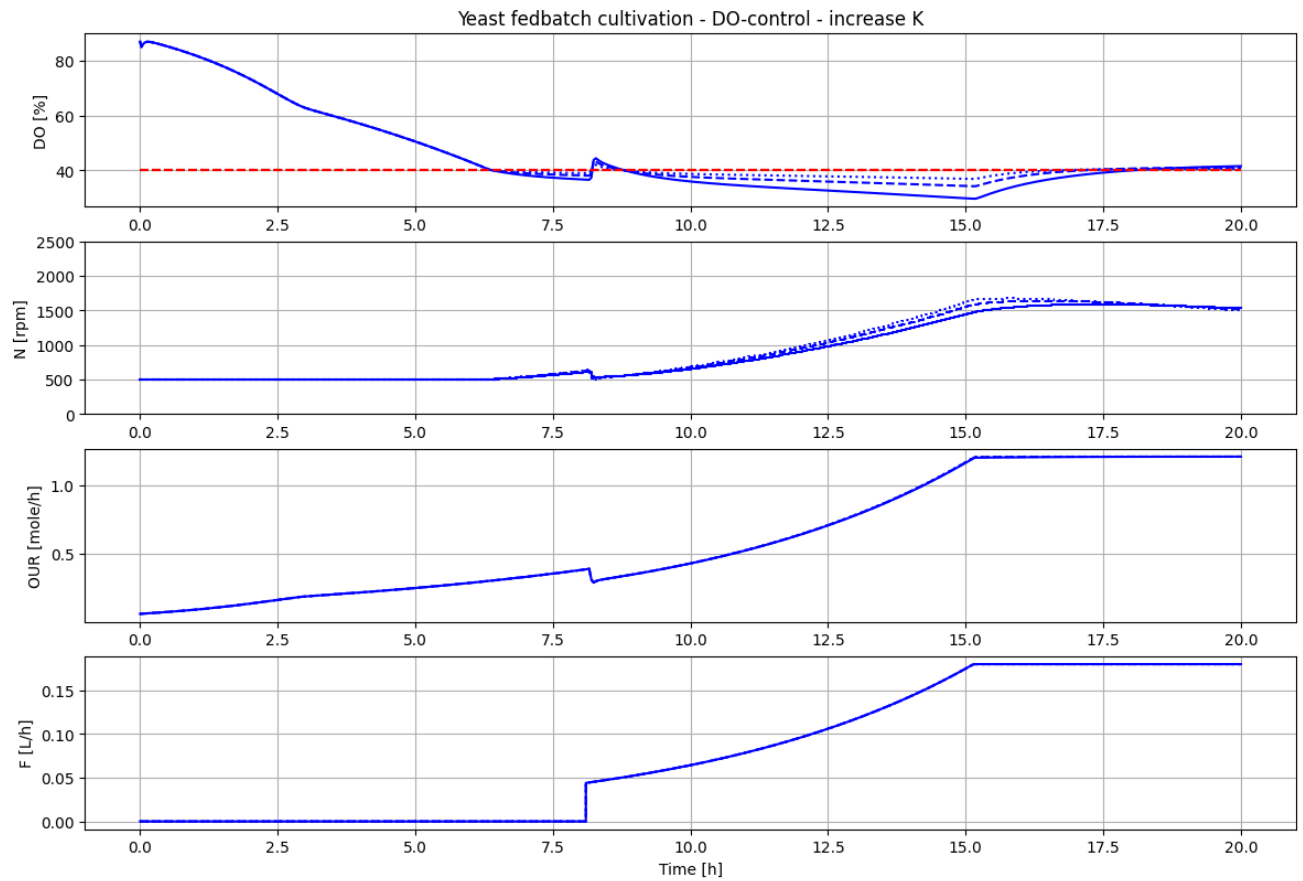


We see that by a higher control gain  $K$  the DO-control error get smaller and the stability of the control system is maintained.

**Exercise** I leave for you to study the impact variation of the  $T_i$ -parameter. Just make a new cell below. Then copy and paste the cell above and change parameter to  $T_i$ .

## ✓ Sensitivity to changes in feed-profile


Now, let us focus on investigating impact of changes in the feed-profile. The goal is to increase the produced cell mass without accumulation of by-product ethanol. Simulation can bring some insight into how behaviour of the different variables change when by-product is formed. This insight can help to interpret experimental results.

```
# Let us check the sensitivity to changes in the feed profile design
newplot(title='Yeast fedbatch cultivation', plotType='Overview')
for value in [0.044, 0.038, 0.050]: par(F_start=value); simu(20)
```

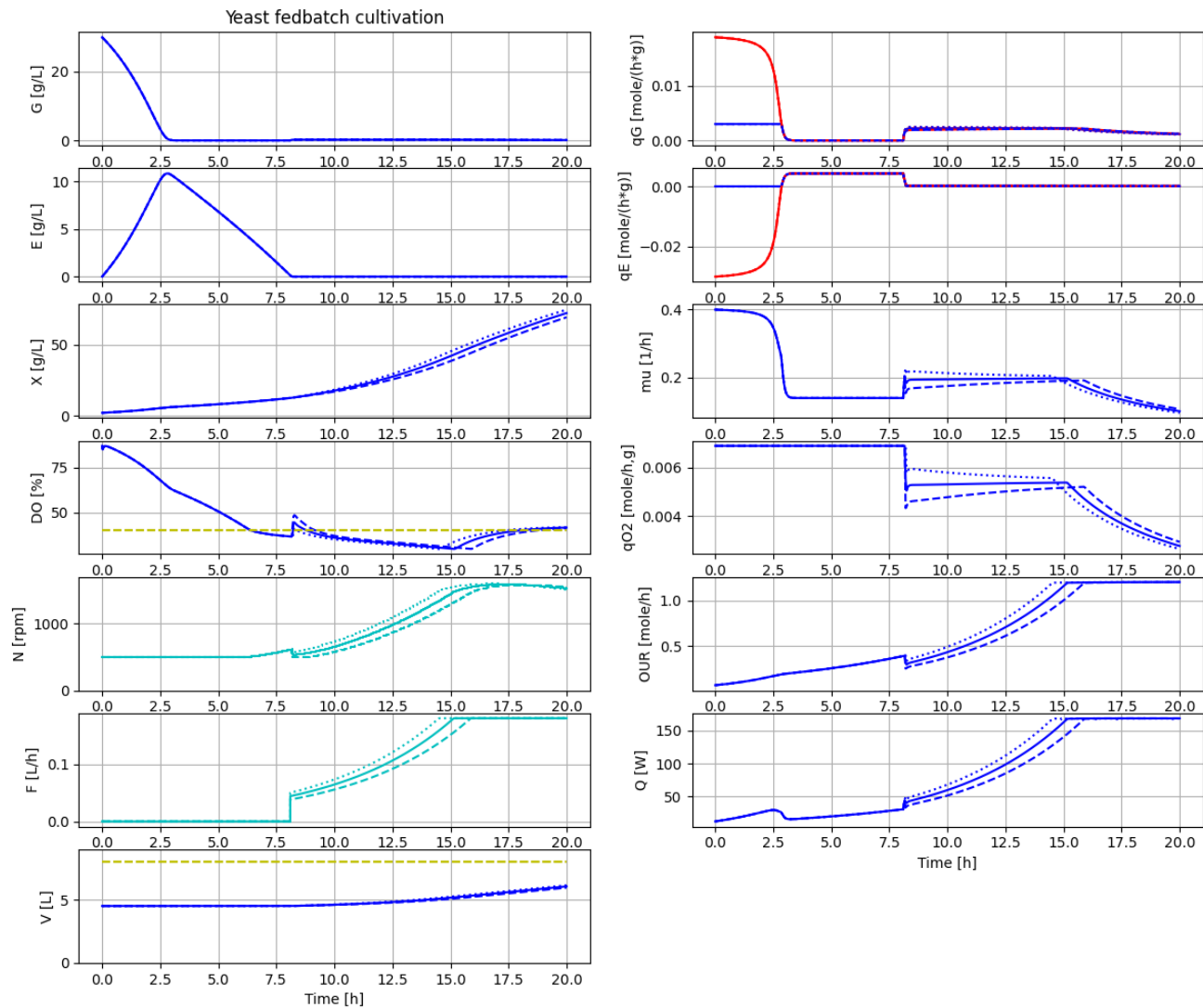


The variation in  $F_{\text{start}}$  has an impact and we see that the actual growth rate during fedbatch phase do converge to the set growth rate of the feed, but it takes more than 5 hours.

```
# Let us investigate a feedprofile that is closer to the maximal capacity
newplot(title='Yeast fedbatch cultivation', plotType='Overview')
par(F_start=0.044, mu_feed=0.20); simu(20)
par(F_start=0.050, mu_feed=0.22); simu(20)
par(F_start=0.057, mu_feed=0.26); simu(20)
```

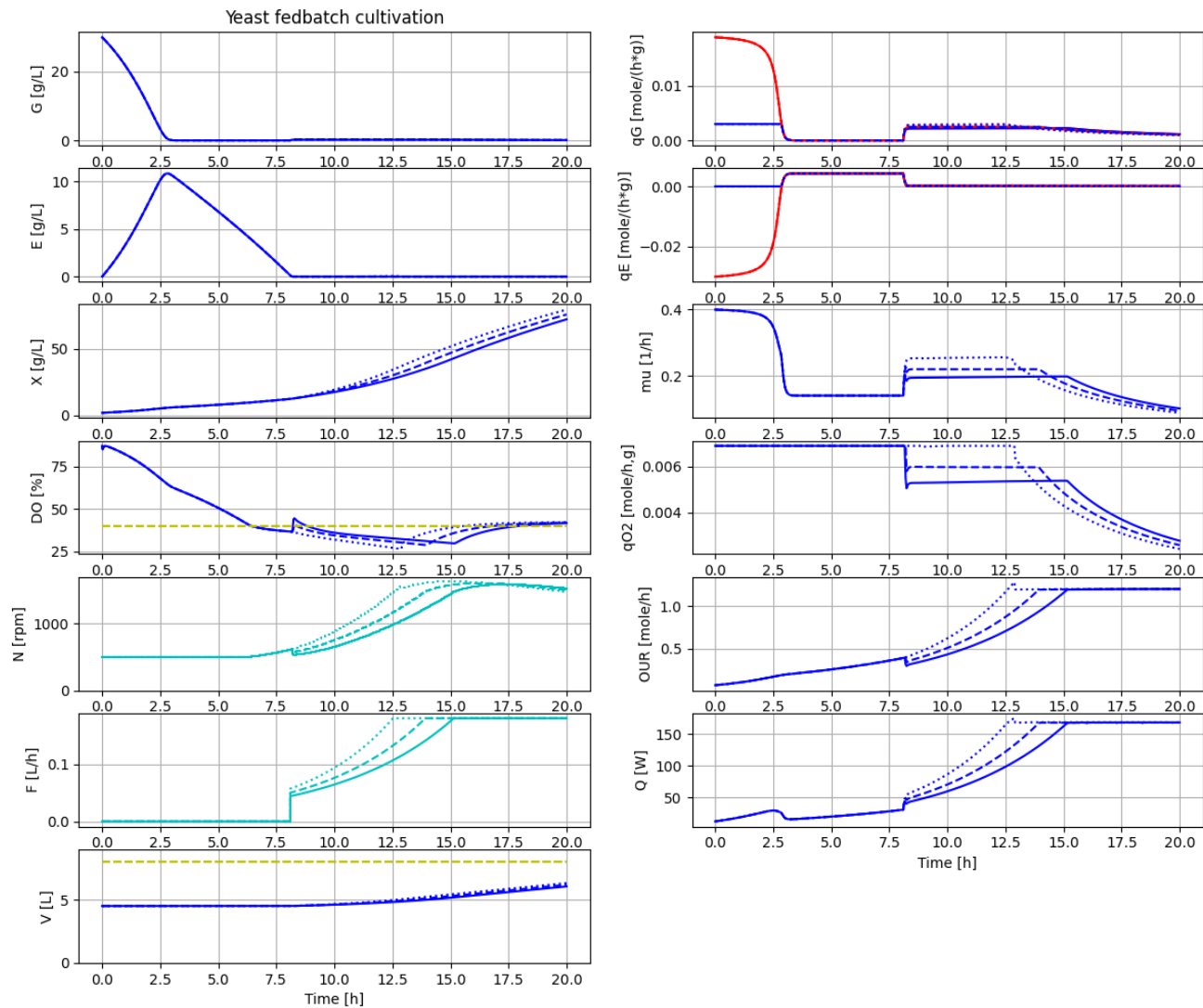

```
# And let us see what happens if the feedprofile exceed the culture capacity  
newplot(title='Yeast fedbatch cultivation', plotType='Overview')  
par(F_start=0.057, mu_feed=0.26); simu(20)  
par(F_start=0.063, mu_feed=0.28); simu(20)  
par(F_start=0.044, mu_feed=0.20)
```



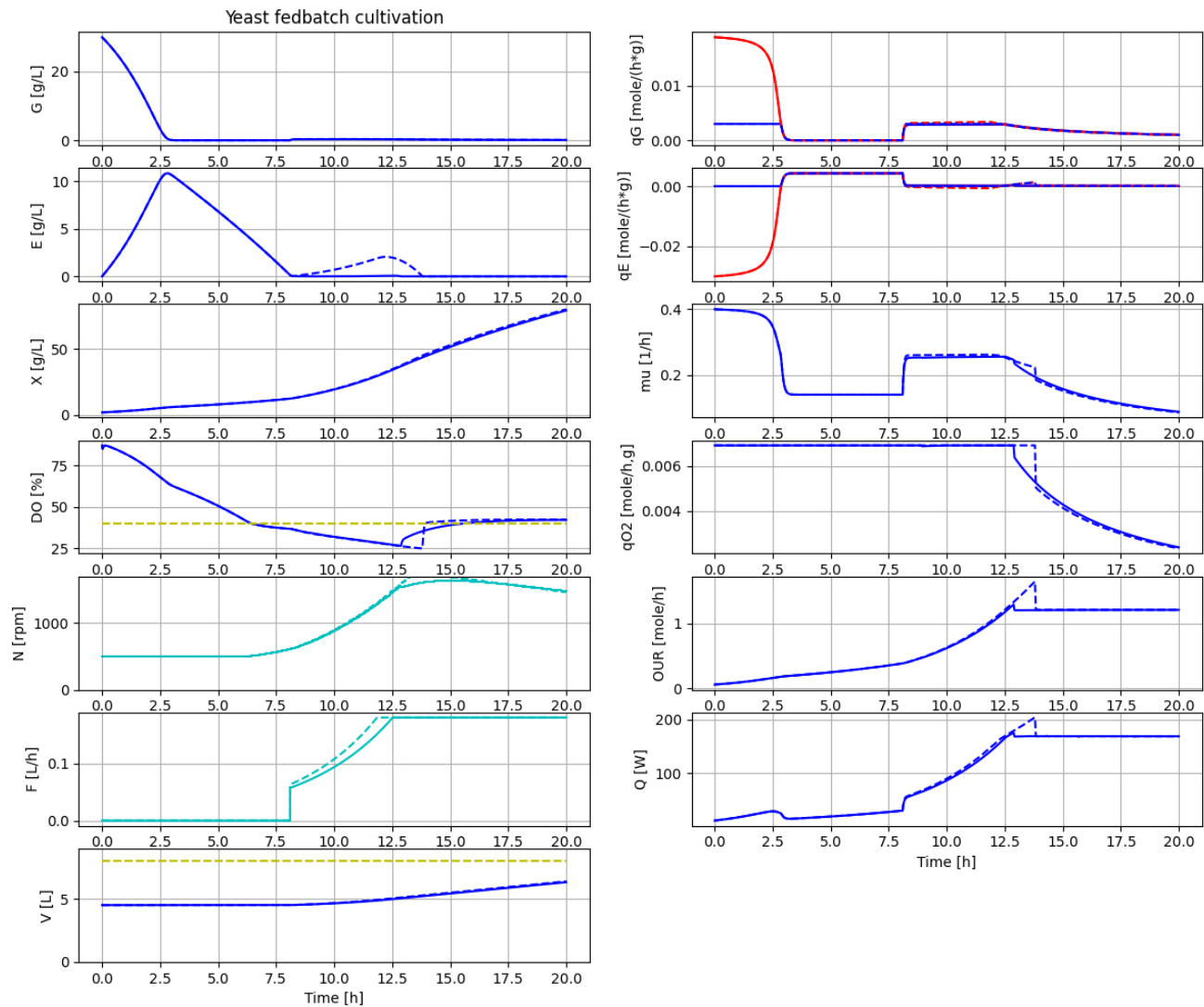

Note that with the feedprofile that exceed culture respiratory capacity, ethanol is accumulated during time 8-12.5 hours. When the feedprofile then is constant from time 12.5 hours and on, then the accumulated ethanol is consumed over about an hour. This leads to a higher oxygen demand and heat production during this time. The specific cell growth rate is also slightly higher during this period.

**Exercise** You can investigate the impact of changing the maximal feedrate  $F_{\max}$ . Make sure that the DO level do not get too low.

```
# Check of simu('cont')
newplot(title='Yeast fedbatch cultivation', plotType='Overview')
par(F_start=0.057, mu_feed=0.26); simu(10)
simu(10, 'cont')
```

